



**ImSpector V8/V10**



**ImSpector V10E**



**ImSpector N17E**

## OPTIONS, FORE OPTICS

- Fore optics, Standard series: OL8, OL12, OL17, OL23 and OL35 for 2/3" or smaller detector.
- Fore optics, Enhanced series: OLEWIDE, OLE18.5, and OLE23 for 2/3" or larger detector. Optimized for Enhanced series.
- Fore optics, OLESMACRO, OLES9, OLES15, OLES22.5, OLES30 and OLES56 for N17E


## OPTIONS, ACCESSORIES

- Mechanical shutter (Enhanced series)
- Collection fiber optics
- Order blocking filters; OBF 570 (rectangular 14 x 12mm or circular 20mm Ø and 17mm Ø) for V10 and V10E


Specim ImSpectors are designed for the VIS (380 - 800 nm), VNIR (400 - 1000 nm) and NIR (900 - 1700 nm) wavelength ranges. These spectrographs provide a straightforward, high performance, yet cost-effective method of integration. When combined with scientific grayscale CCD or CMOS cameras or InGaAs sensor, the combination provides a line-scan Spectral Imaging device.

## DIMENSIONS


### ImSpector V8/V10

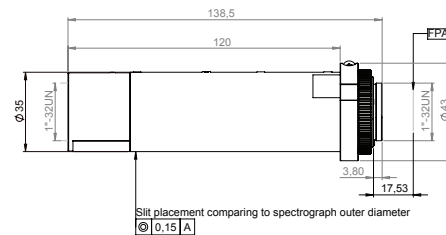

### ImSpector V10E

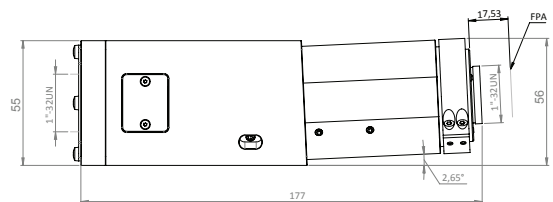

### ImSpector N17E

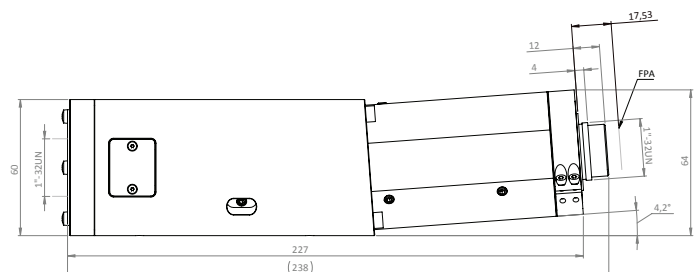

ImSpector	V8	V10	V10E	N17E
Optical characteristics				
Spectral range	380 - 800 nm *1	400 - 1000 nm *1	400 - 1000 nm *1	900 - 1700 nm
Dispersion	66 nm / mm	93.9 nm / mm	97.5 nm / mm	110 nm / mm
Spectral resolution	6 nm (with 80 µm slit)	9 nm (with 80 µm slit)	2.8 nm (with 30 µm slit)	5 nm (with 30 µm slit)
Image size	6.4 (spectral) x 8.8 (spatial) mm corresponding to standard ⅔" image sensor	6.4 (spectral) x 8.8 (spatial) mm corresponding to standard ⅔" image sensor	Max 6.15 (spectral) x 14.2 (spatial) mm	Max 7.6 (spectral) x 14.2 (spatial) mm
Spatial resolution	Rms spot radius < 30 µm	Rms spot radius < 30 µm	Rms spot radius < 9 µm	Rms spot radius < 15 µm
Aberrations	Insignificant astigmatism	Insignificant astigmatism	No astigmatism	No astigmatism
Bending of spectral lines across spatial axis	Smile < 45 µm	Smile < 45 µm	Smile < 1.5 µm	Smile < 5 µm
Bending of spatial lines across spectral axis	Keystone < 20 µm	Keystone < 20 µm	Keystone < 1 µm	Keystone < 5 µm
Working F-number	F/2.8	F/2.8	F/2.4	F/2.0
Slit width, default	50 µm (30 and 80 µm on request)	50 µm (30 and 80 µm on request)	30 µm (50 and 80 µm on request)	30 µm
Slit length	8.8 mm	8.8 mm	14.2 mm	14.2 mm
Optical input	N/A	N/A	Telecentric	Telecentric
Efficiency	> 50% independent of polarization			
Mechanical characteristics				
Size	D 35 x L 139 mm	D 35 x L 139 mm	W 60 x H 60 x L 175 mm	W 60 x H 60 x L 220 mm
Weight	300 g	300 g	1100 g	1500 g
Body	Anonized aluminium tube			
Lens and camera mount	Standard C-mount adapter			
User adjustments	Image axis relative to detector rows, adjustable back focal length +/- 1mm			
Environmental characteristics				
Storage	-20 ... +85 °C			
Operating	+5 ... +40 °C, non-condensing			

\*1 Order blocking filter is available for mounting in front of the detector window.

Specim, Spectral Imaging Ltd. ▪ A Konica Minolta Company ▪ POB 110, FI-90591 Oulu Finland ▪ Elektroniikkatie 13, Oulu Finland  
Tel +358 (0) 10 4244 400 ▪ VAT Identification number FI10079234 ▪ info@specim.com ▪ www.specim.com

All specifications provided are typical values and are subject to change without notice. These values are not guaranteed and should not be considered absolute limits. Specim, Spectral Imaging Ltd. reserves the right to change or improve its products and specifications and to make changes in content without obligation to notify any person or organization of such changes or improvements.